

Heated Wire Circuits



YOUR BEST PARTNER



DA CHUNG MEDICAL CO., LTD.

Manufacturer site: No.600, Sec.1, Dong Shan Rd., Dong Shan Township, Yilan County, Taiwan
TEL: +886 3 959 4777 FAX: +886 3 959 4707
Taipei office: Rm. 1420, 14F., No. 206, Sec. 1, Keelung Rd., Xinyi Dist., Taipei City, Taiwan
TEL: +886 2 6600 4520 E-mail: info@dc-med.com.tw




D.C. Medical heated wire circuit as an excellent solution to respiratory humidification. Full series compatible with F&P MR250 (and warmers with similar design). Unique and dedicated design ensure ease of use for healthcare professionals, and mean time, secure safety and efficacy for patients



Unique water trap design

- Circuit auto-seal mechanism, care-free water pouring
- Less ventilator alarm^a
- Less opportunistic infection^b


a. Compared with non auto-seal water trap model
 b. Respiratory Care May 2019, 64 (5) 545-549

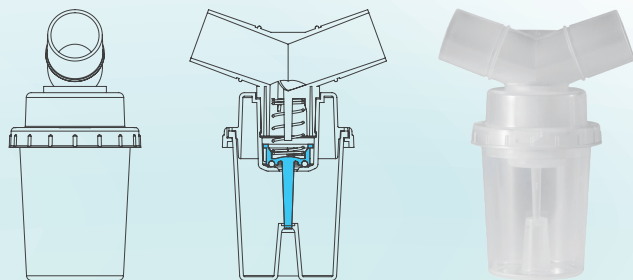

Patented auto-refill chamber

- Automatically keep water in a steady level^a
- 100% tested during production, hassle free, load and forget^a
- PP materil, zero BPA, zero interference to human body, safe for all ages


a. DC Medical, production and inspection standard operation process

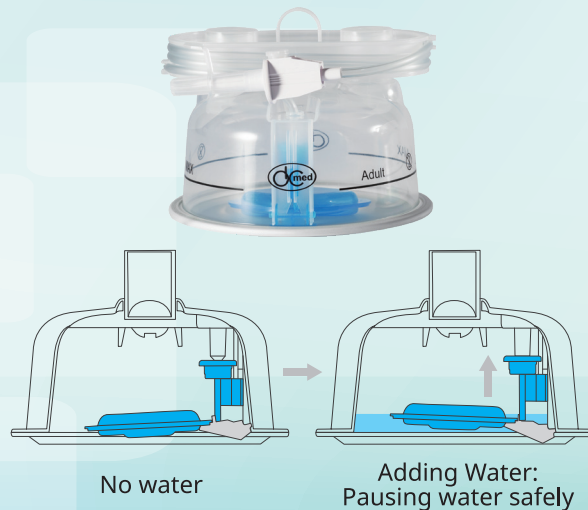

Optimal wire length

- Efficient and stable heating
- Superior heating function, 21% faster in reaching 37 degree celcius^a
- Less condensation inside circuit, 20% reduction in the volume of water^b


a. Compared with economic alternative
 b. Respiratory Care May 2019, 64 (5) 545-549

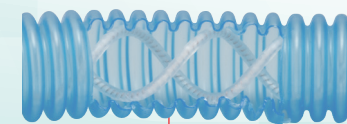

Unique combination of material


- 21%-44% less expansion against pressure^a
- Air-tight without adhesive/sealant, zero leakage under 60mmH2O^b

a. Compared with common PP circuit
 b. DC Medical product master file

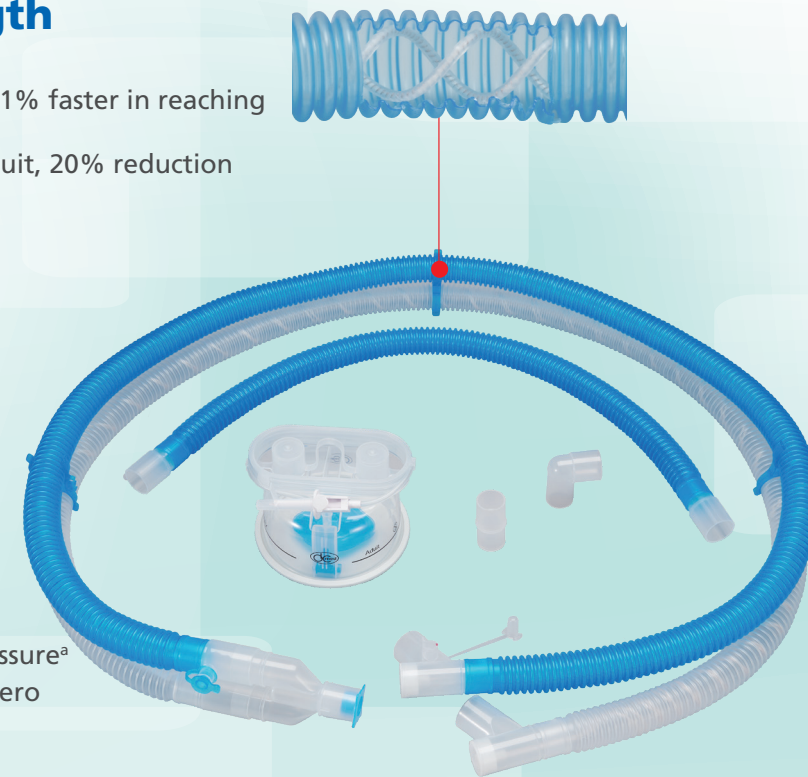


ADULT

- ① Heated expiratory limb: 150cm
- ② Non-heated expiratory limb: 90 cm + 60cm
- ③ Inspiratory limb: 150 cm, 22ID, MR850 system chamber adaptor
- ④ Chamber limb: 60 cm
- ⑤ 55ml Water trap 22 mmOD
- ⑥ Auto-fill chamber connection ports: 22 mmOD
- ⑦ Connector kit: straight connector, 22mmOD; elbow connector, 22 mmOD, 22 mmID
- ⑧ Temperature port: 7.6 mmID
- ⑨ Patient wye connector :22/15 mm

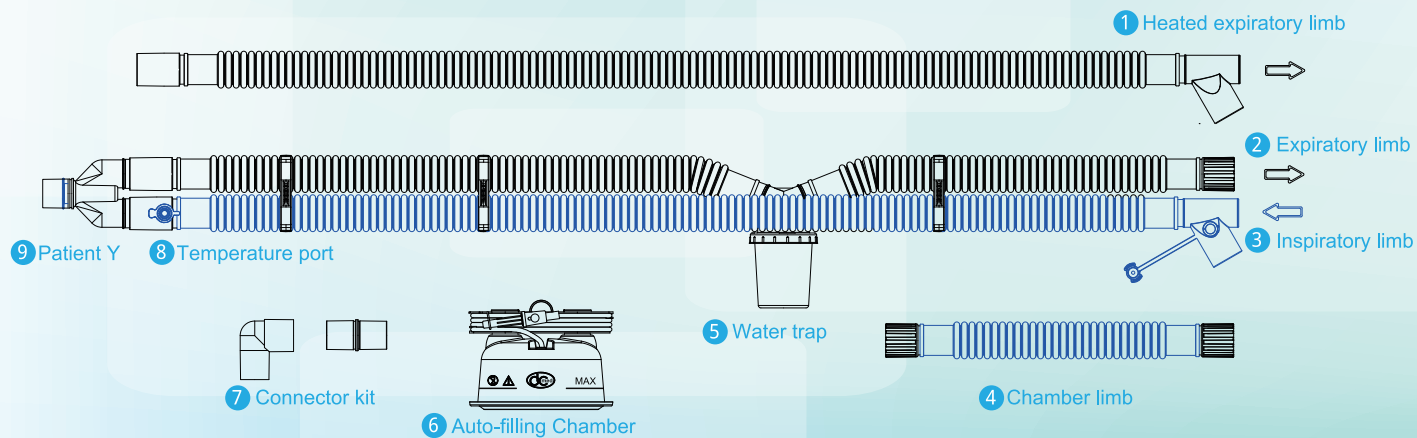

NEONATE

- ① Heated expiratory limb: 140 cm
- ② Non-heated expiratory limb: 90+60 cm
- ③ Inspiratory limb: 110+30 cm, 22 mmID, MR 850 style humidifier adaptor
- ④ Chamber limb: 60 cm
- ⑤ 55ml Water trap 10 mmOD
- ⑥ Auto-fill chamber connection ports: 22 mmOD
- ⑦ Connector kit
- ⑧ Temperature port: 7.6 mmID
- ⑨ Patient wye connector: 15 mmID
- ⑩ Pressure port

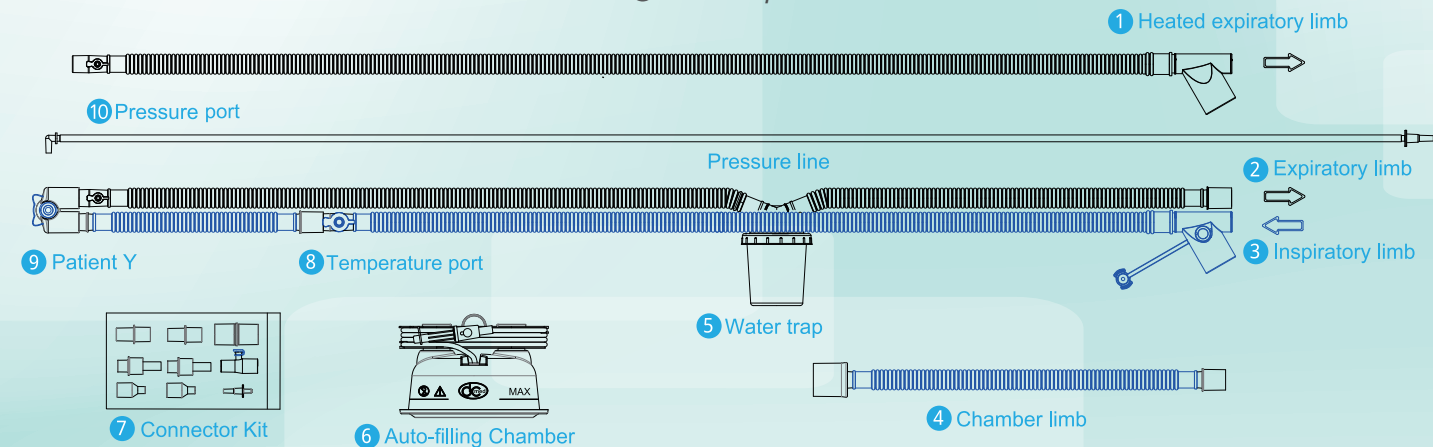

ADULT

REF	Heated Wire	Description	Package	Accessory
3410	Insp.	Single Heated Wire Circuit	1 set/bag, 10 sets/carton	Chamber included
3420	Insp. & Exp.	Dual Heated Wire Circuit	1 set/bag, 10 sets/carton	Chamber included
3411	Insp.	Single Heated Wire Circuit	1 set/bag, 10 sets/carton	
3421	Insp. & Exp.	Dual Heated Wire Circuit	1 set/bag, 10 sets/carton	
3413	Insp.	Inspiratory Heated Wire Circuit	1 set/bag, 15 sets/carton	Chamber included
3490		Auto-filling Humidification Chamber	1 set/bag, 40 sets/carton	
3493		Manual-filling Humidification Chamber	1 set/bag, 40 sets/carton	

NEONATE

REF	Heated Wire	Description	Package	Accessory
3440	Insp.	Single Heated Wire Circuit	1 set/bag, 20 sets/carton	Chamber, Connector Kit
3450	Insp. & Exp.	Dual Heated Wire Circuit	1 set/bag, 20 sets/carton	Chamber, Connector Kit
3441	Insp.	Single Heated Wire Circuit	1 set/bag, 20 sets/carton	Connector Kit
3451	Insp. & Exp.	Dual Heated Wire Circuit	1 set/bag, 20 sets/carton	Connector Kit
3445	Insp.	Single-heated Wire Circuit with Chamber (Auto-filling) for SLE5000	1 set/bag, 20 sets/carton	Chamber, Connector Kit
3449	Insp.	Single-heated Wire Circuit with Chamber (Manual-filling) for SLE5000	1 set/bag, 20 sets/carton	Chamber, Connector Kit
3446	Insp.	Single Heated Wire Circuit for SLE5000	1 set/bag, 20 sets/carton	Connector Kit
3443	Insp.	Inspiratory Heated Wire Circuit with Chamber (Auto-filling)	1 set/bag, 20 sets/carton	Chamber, Connector Kit
3484		SLE5000 Connector (Restrictor)	1000 pcs/bag, 1 bag/carton	Bulk package
3495		Auto-filling Humidification Chamber	1 set/bag, 40 sets/carton	
3498		Manual-filling Humidification Chamber	1 set/bag, 40 sets/carton	